

Unnamed Protocol

OT-2 Simulation Protocol Report

Summary

Pipetting steps:	138
Labware items:	3
Tracked plates:	2
Deck snapshots:	139


Labware Layout













Slot 1:	Opentrons OT-2 96 Tip Rack 30...
Slot 2:	NEST 96 Well Plate 200 μ L Flat
Slot 3:	NEST 96 Well Plate 200 μ L Flat

Deck Layout

Slot 10 -	Slot 11 -	Trash -
Slot 7 -	Slot 8 -	Slot 9 -
Slot 4 -	Slot 5 -	Slot 6 -
Slot 1 Opentrons OT-2 96 Tip R...	Slot 2 NEST 96 Well Plate 200 ...	Slot 3 NEST 96 Well Plate 200 ...

Unnamed Protocol · Initial State

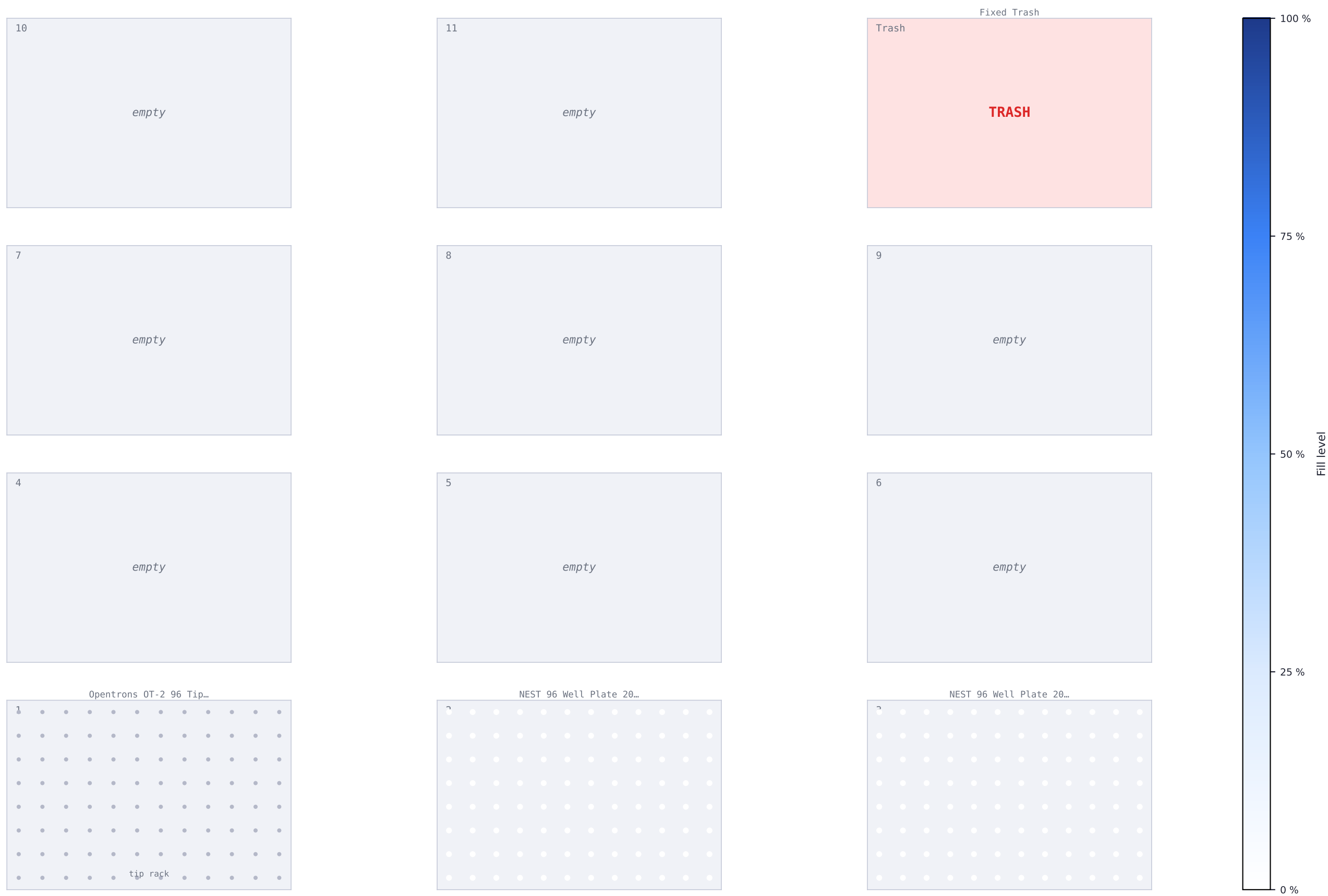

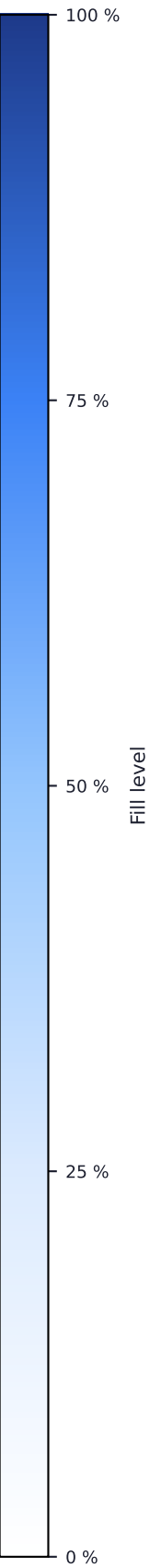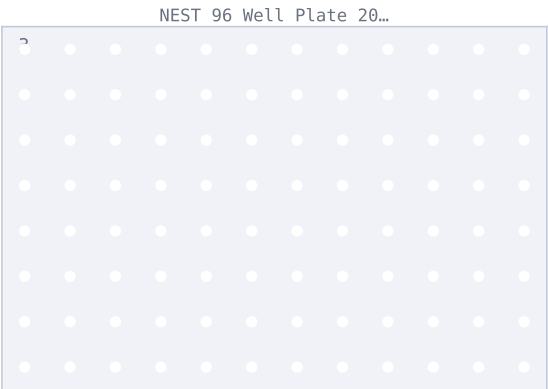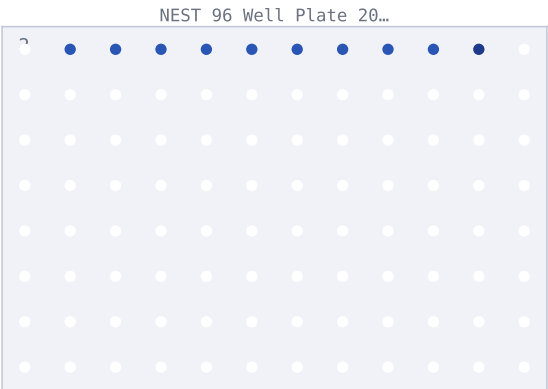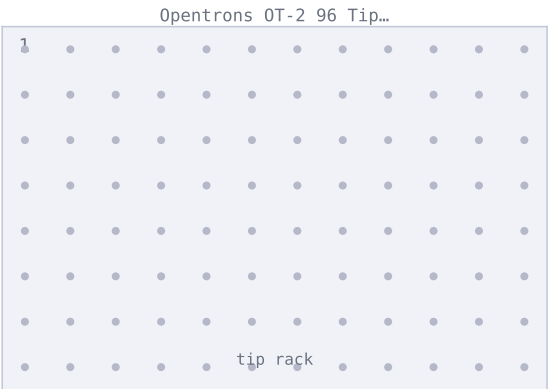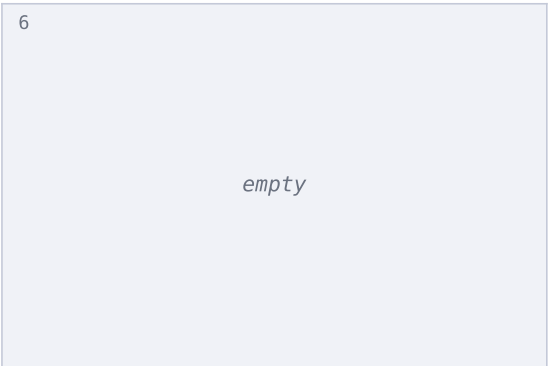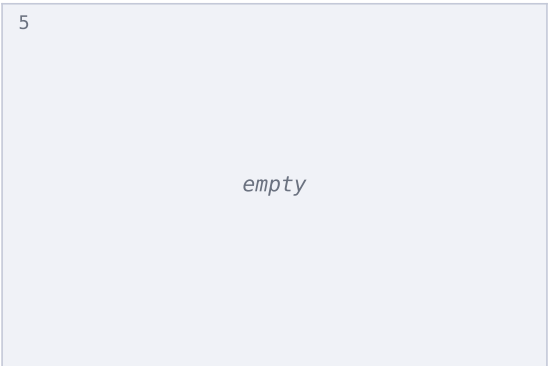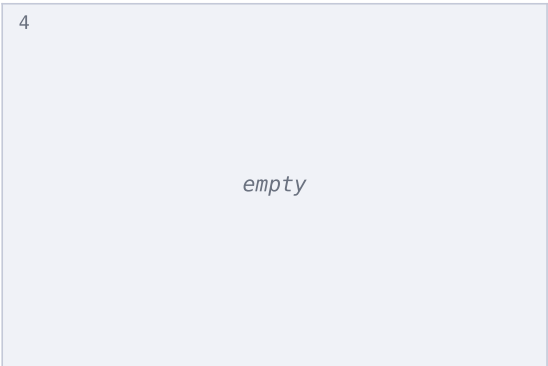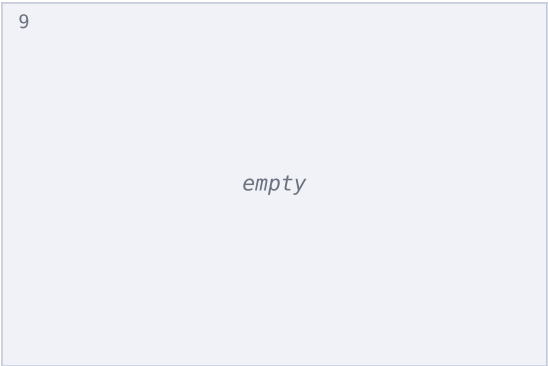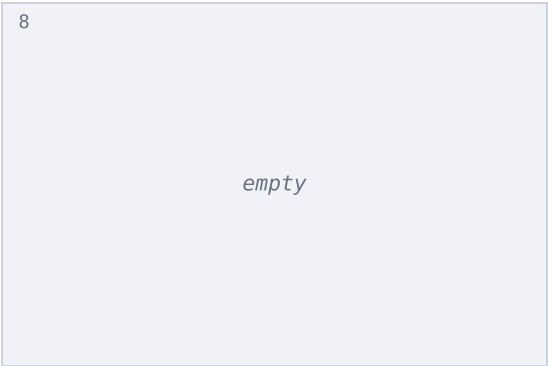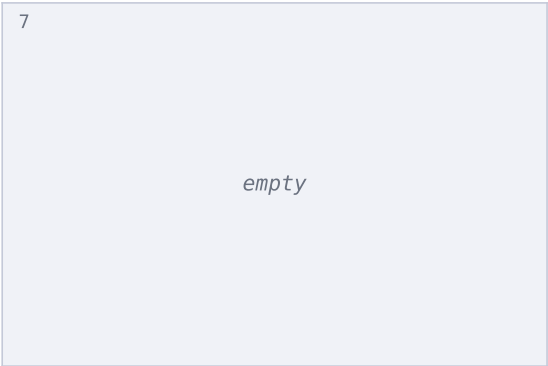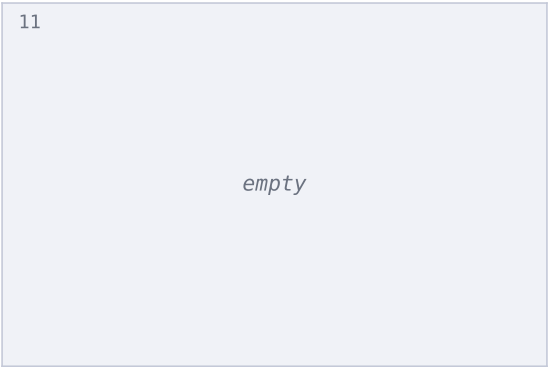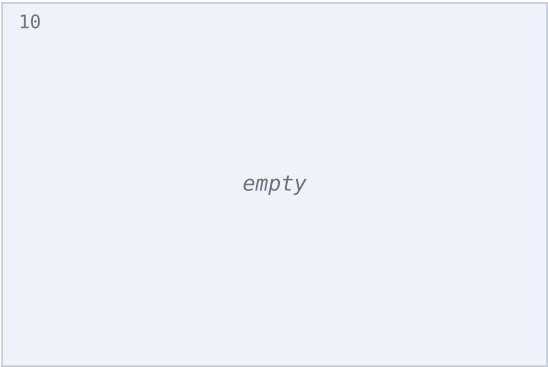

Initial StateFinal State

Plate Comparison – Slot 2: NEST 96 Well Plate 200 μ L Flat

1 2 3 4 5 6 7 8 9 10 11 12

A

B

C

D

E

F

G

H

1 2 3 4 5 6 7 8 9 10 11 12

A

B

C

D

E

F

G

H

180

180

180

180

180

180

180

180

180

200

200

150

100

50

0


Volume (μ L) / max = 200 μ L

Initial StateFinal State

Plate Comparison – Slot 3: NEST 96 Well Plate 200 μ L Flat



1 2 3 4 5 6 7 8 9 10 11 12

A
B
C
D
E
F
G
H



1 2 3 4 5 6 7 8 9 10 11 12

A
B
C
D
E
F
G
H

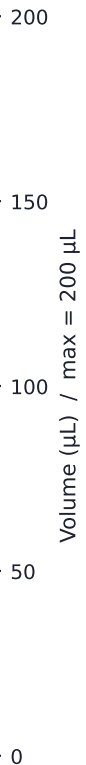

Step Log (page 1 of 4)

Unnamed Protocol

Step	Action	Vol (μL)	Slot	Well	Labware
0	Pick Up Tip	—	1	A1	Opentrons OT-2 96 Tip Rack 300 μL
1	Note	—	—	—	—
2	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
3	Dispense	180.0	2	A2	NEST 96 Well Plate 200 μL Flat
4	Note	—	—	—	—
5	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
6	Dispense	180.0	2	A3	NEST 96 Well Plate 200 μL Flat
7	Note	—	—	—	—
8	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
9	Dispense	180.0	2	A4	NEST 96 Well Plate 200 μL Flat
10	Note	—	—	—	—
11	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
12	Dispense	180.0	2	A5	NEST 96 Well Plate 200 μL Flat
13	Note	—	—	—	—
14	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
15	Dispense	180.0	2	A6	NEST 96 Well Plate 200 μL Flat
16	Note	—	—	—	—
17	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
18	Dispense	180.0	2	A7	NEST 96 Well Plate 200 μL Flat
19	Note	—	—	—	—
20	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
21	Dispense	180.0	2	A8	NEST 96 Well Plate 200 μL Flat
22	Note	—	—	—	—
23	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
24	Dispense	180.0	2	A9	NEST 96 Well Plate 200 μL Flat
25	Note	—	—	—	—
26	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
27	Dispense	180.0	2	A10	NEST 96 Well Plate 200 μL Flat
28	Note	—	—	—	—
29	Aspirate	180.0	3	A1	NEST 96 Well Plate 200 μL Flat
30	Dispense	180.0	2	A11	NEST 96 Well Plate 200 μL Flat
31	Drop Tip	—	12	—	A1 of Opentrons Fixed Trash
32	Pick Up Tip	—	1	B1	Opentrons OT-2 96 Tip Rack 300 μL
33	Note	—	—	—	—
34	Aspirate	20.0	2	A1	NEST 96 Well Plate 200 μL Flat

Step Log (page 2 of 4)

Unnamed Protocol

Step	Action	Vol (μL)	Slot	Well	Labware
35	Dispense	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
36	Mix	—	—	—	—
37	Aspirate	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
38	Dispense	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
39	Aspirate	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
40	Dispense	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
41	Aspirate	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
42	Dispense	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
43	Note	—	—	—	—
44	Aspirate	20.0	2	A2	NEST 96 Well Plate 200 μL Flat
45	Dispense	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
46	Mix	—	—	—	—
47	Aspirate	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
48	Dispense	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
49	Aspirate	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
50	Dispense	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
51	Aspirate	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
52	Dispense	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
53	Note	—	—	—	—
54	Aspirate	20.0	2	A3	NEST 96 Well Plate 200 μL Flat
55	Dispense	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
56	Mix	—	—	—	—
57	Aspirate	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
58	Dispense	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
59	Aspirate	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
60	Dispense	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
61	Aspirate	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
62	Dispense	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
63	Note	—	—	—	—
64	Aspirate	20.0	2	A4	NEST 96 Well Plate 200 μL Flat
65	Dispense	20.0	2	A5	NEST 96 Well Plate 200 μL Flat
66	Mix	—	—	—	—
67	Aspirate	20.0	2	A5	NEST 96 Well Plate 200 μL Flat
68	Dispense	20.0	2	A5	NEST 96 Well Plate 200 μL Flat
69	Aspirate	20.0	2	A5	NEST 96 Well Plate 200 μL Flat

Step Log (page 3 of 4)

Unnamed Protocol

Step	Action	Vol (μL)	Slot	Well	Labware
70	Dispense	20.0	2	A5	NEST 96 Well Plate 200 μL Flat
71	Aspirate	20.0	2	A5	NEST 96 Well Plate 200 μL Flat
72	Dispense	20.0	2	A5	NEST 96 Well Plate 200 μL Flat
73	Note	—	—	—	—
74	Aspirate	20.0	2	A5	NEST 96 Well Plate 200 μL Flat
75	Dispense	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
76	Mix	—	—	—	—
77	Aspirate	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
78	Dispense	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
79	Aspirate	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
80	Dispense	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
81	Aspirate	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
82	Dispense	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
83	Note	—	—	—	—
84	Aspirate	20.0	2	A6	NEST 96 Well Plate 200 μL Flat
85	Dispense	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
86	Mix	—	—	—	—
87	Aspirate	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
88	Dispense	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
89	Aspirate	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
90	Dispense	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
91	Aspirate	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
92	Dispense	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
93	Note	—	—	—	—
94	Aspirate	20.0	2	A7	NEST 96 Well Plate 200 μL Flat
95	Dispense	20.0	2	A8	NEST 96 Well Plate 200 μL Flat
96	Mix	—	—	—	—
97	Aspirate	20.0	2	A8	NEST 96 Well Plate 200 μL Flat
98	Dispense	20.0	2	A8	NEST 96 Well Plate 200 μL Flat
99	Aspirate	20.0	2	A8	NEST 96 Well Plate 200 μL Flat
100	Dispense	20.0	2	A8	NEST 96 Well Plate 200 μL Flat
101	Aspirate	20.0	2	A8	NEST 96 Well Plate 200 μL Flat
102	Dispense	20.0	2	A8	NEST 96 Well Plate 200 μL Flat
103	Note	—	—	—	—
104	Aspirate	20.0	2	A8	NEST 96 Well Plate 200 μL Flat

Step Log (page 4 of 4)

Unnamed Protocol

Step	Action	Vol (μL)	Slot	Well	Labware
105	Dispense	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
106	Mix	—	—	—	—
107	Aspirate	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
108	Dispense	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
109	Aspirate	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
110	Dispense	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
111	Aspirate	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
112	Dispense	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
113	Note	—	—	—	—
114	Aspirate	20.0	2	A9	NEST 96 Well Plate 200 μL Flat
115	Dispense	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
116	Mix	—	—	—	—
117	Aspirate	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
118	Dispense	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
119	Aspirate	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
120	Dispense	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
121	Aspirate	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
122	Dispense	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
123	Note	—	—	—	—
124	Aspirate	20.0	2	A10	NEST 96 Well Plate 200 μL Flat
125	Dispense	20.0	2	A11	NEST 96 Well Plate 200 μL Flat
126	Mix	—	—	—	—
127	Aspirate	20.0	2	A11	NEST 96 Well Plate 200 μL Flat
128	Dispense	20.0	2	A11	NEST 96 Well Plate 200 μL Flat
129	Aspirate	20.0	2	A11	NEST 96 Well Plate 200 μL Flat
130	Dispense	20.0	2	A11	NEST 96 Well Plate 200 μL Flat
131	Aspirate	20.0	2	A11	NEST 96 Well Plate 200 μL Flat
132	Dispense	20.0	2	A11	NEST 96 Well Plate 200 μL Flat
133	Drop Tip	—	12	—	A1 of Opentrons Fixed Trash
134	Note	—	—	—	—
135	Note	—	—	—	—
136	Note	—	—	—	—
137	Note	—	—	—	—